


For short implant placement, implant major diameter needs to be \leq the bur (average diameter) at the 8mm laser mark. Please refer to page 16 in the Instructions for Use Manual.

In Ridge Expansion cases, please oversize your osteotomy and make sure that the crest diameter is equal to or larger than the implant major diameter.

In Hard Bone (Mandible), after Finishing the Full Osteotomy Preparation, Use the Next Larger Size Densah Bur to the 3mm Laser-Mark Depth to make sure the Osteotomy Crestal Diameter is Equal to or Larger than the Implant Major (Crestal) Diameter.

Use Densah Burs in full-step increments for Sinus Lift cases. Example: 2.0mm, 3.0mm, 4.0mm, 5.0mm

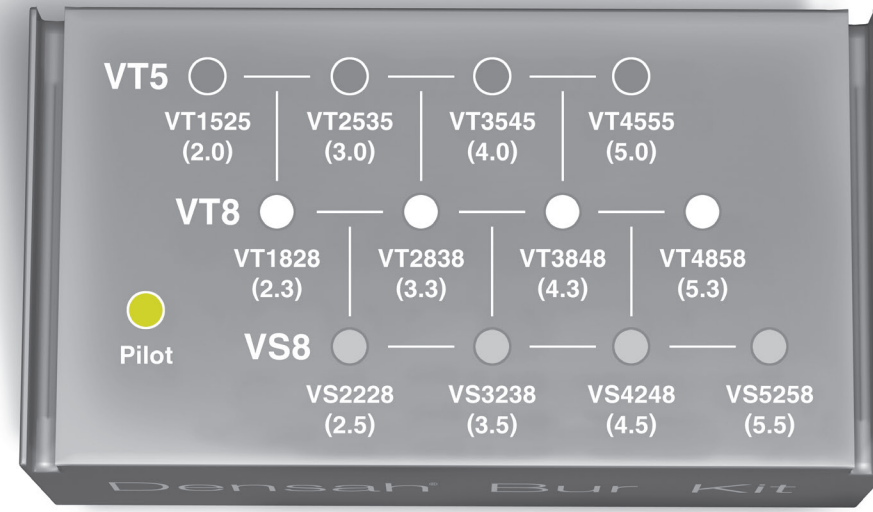

Use large block display to compare Bur identification system when using the schematic below for proper Bur usage

● VT5 Set ○ VT8 Set ● VS8 Set

Densifying Mode CCW (800-1500) RPMs / Cutting Mode CW (800-1500) RPMs

Tag			Axis															
			Soft Bone							Hard Bone (Mandible)								
			In densifying mode make sure your osteotomy is 1.0 mm deeper than the actual implant final length. In extreme hard bone, utilize DAC (Densify After Cut) Protocol. Find protocol in IFU.															
Geometry	Major Ø	Minor Ø	Pilot	Bur 1	Bur 2	Bur 3	Bur 4	Bur 5	Densah® Bur Block Display	Pilot	Bur 1	Bur 2	Bur 3	Bur 4	Bur 5	Bur 6	Bur 7	Densah® Bur Block Display
Straight	3.3		Pilot	VT1828 (2.3)	VS2228* (2.5)	—	—	—		Pilot	VT1525 (2.0)	VT1828 (2.3)	VS2228* (2.5)	—	—	—	—	
Straight	3.75		Pilot	VT1828 (2.3)	VT2838 (3.3)	VS3238* (3.5)	—	—		Pilot	VT1828 (2.3)	VT2838 (3.3)	VS3238* (3.5)	—	—	—	—	
Straight	4.2		Pilot	VT1828 (2.3)	VT2838 (3.3)	VS3238* (3.5)	—	—		Pilot	VT1828 (2.3)	VT2838 (3.3)	VS3238* (3.5)	—	—	—	—	
Straight	5.0		Pilot	VT1828 (2.3)	VT2838 (3.3)	VT3848 (4.3)	VS4248* (4.5)	—		Pilot	VT1828 (2.3)	VT2535 (3.0)	VT2838 (3.3)	VT3545 (4.0)	VT3848 (4.3)	VS4248* (4.5)	—	
Straight	6.0		Pilot	VT1828 (2.3)	VT2838 (3.3)	VT3848 (4.3)	VT4858 (5.3)	VS5258* (5.5)		Pilot	VT1828 (2.3)	VT2838 (3.3)	VT3545 (4.0)	VT3848 (4.3)	VT4555 (5.0)	VT4858 (5.3)	VS5258* (5.5)	

*Denotes implant placement.

*Clinician experience and judgment should be used in conjunction with the Densifying Reference Guide recommendation and suggested use protocols.

For short implant placement, implant major diameter needs to be \leq the bur (average diameter) at the 8mm laser mark. Please refer to page 16 in the Instructions for Use Manual.

In Ridge Expansion cases, please oversize your osteotomy and make sure that the crest diameter is equal to or larger than the implant major diameter.

In Hard Bone (Mandible), after Finishing the Full Osteotomy Preparation, Use the Next Larger Size Densah Bur to the 3mm Laser-Mark Depth to make sure the Osteotomy Crestal Diameter is Equal to or Larger than the Implant Major (Crestal) Diameter.

Use Densah Burs in full-step increments for Sinus Lift cases. Example: 2.0mm, 3.0mm, 4.0mm, 5.0mm



Use large block display to compare Bur identification system when using the schematic below for proper Bur usage

● VT5 Set ○ VT8 Set ● VS8 Set

Densifying Mode CCW (800-1500) RPMs / Cutting Mode CW (800-1500) RPMs

Tag		Massif																	
		Soft Bone								Hard Bone (Mandible)									
Geometry	Major Ø	Minor Ø	Pilot	Bur 1	Bur 2	Bur 3	Bur 4	Bur 5	Densah® Bur Block Display	Pilot	Bur 1	Bur 2	Bur 3	Bur 4	Bur 5	Bur 6	Bur 7	Densah® Bur Block Display	
										In densifying mode make sure your osteotomy is 1.0 mm deeper than the actual implant final length. In extreme hard bone, utilize DAC (Densify After Cut) Protocol. Find protocol in IFU.									
Straight	3.75		Pilot	VT1828 (2.3)	VT2838 (3.3)	VS3238* (3.5)	—	—		Pilot	VT1828 (2.3)	VT2838 (3.3)	VS3238* (3.5)	—	—	—	—		
Straight	4.2		Pilot	VT1828 (2.3)	VT2838 (3.3)	VS3238* (3.5)	—	—		Pilot	VT1828 (2.3)	VT2838 (3.3)	VS3238* (3.5)	—	—	—	—		
Straight	5.0		Pilot	VT1828 (2.3)	VT2838 (3.3)	VT3848 (4.3)	VS4248* (4.5)	—		Pilot	VT1828 (2.3)	VT2535 (3.0)	VT2838 (3.3)	VT3545 (4.0)	VT3848 (4.3)	VS4248* (4.5)	—		
Straight	6.0		Pilot	VT1828 (2.3)	VT2838 (3.3)	VT3848 (4.3)	VT4858 (5.3)	VS5258* (5.5)		Pilot	VT1828 (2.3)	VT2838 (3.3)	VT3545 (4.0)	VT3848 (4.3)	VT4555 (5.0)	VT4858 (5.3)	VS5258* (5.5)		

*Denotes implant placement.


*Clinician experience and judgment should be used in conjunction with the Densifying Reference Guide recommendation and suggested use protocols.

For short implant placement, implant major diameter needs to be \leq the bur (average diameter) at the 8mm laser mark. Please refer to page 16 in the Instructions for Use Manual.

In Ridge Expansion cases, please oversize your osteotomy and make sure that the crest diameter is equal to or larger than the implant major diameter.

In Hard Bone (Mandible), after Finishing the Full Osteotomy Preparation, Use the Next Larger Size Densah Bur to the 3mm Laser-Mark Depth to make sure the Osteotomy Crestal Diameter is Equal to or Larger than the Implant Major (Crestal) Diameter.

Use Densah Burs in full-step increments for Sinus Lift cases. Example: 2.0mm, 3.0mm, 4.0mm, 5.0mm

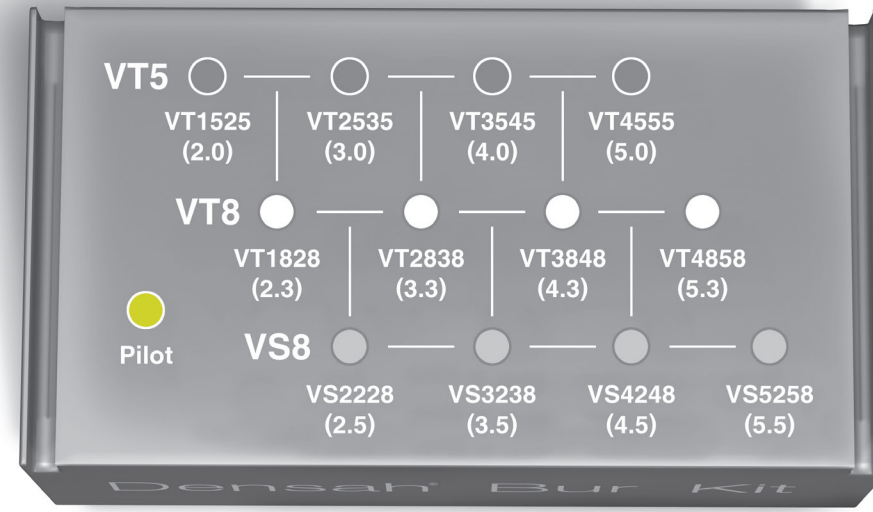

Use large block display to compare Bur identification system when using the schematic below for proper Bur usage

● VT5 Set ○ VT8 Set ● VS8 Set

Densifying Mode CCW (800-1500) RPMs / Cutting Mode CW (800-1500) RPMs

Tag			Crestone														
			Soft Bone						Hard Bone (Mandible)								
Geometry	Major Ø	Minor Ø	Pilot	Bur 1	Bur 2	Bur 3	Bur 4	Densah® Bur Block Display	In densifying mode make sure your osteotomy is 1.0 mm deeper than the actual implant final length. In extreme hard bone, utilize DAC (Densify After Cut) Protocol. Find protocol in IFU.								Densah® Bur Block Display
									Pilot	Bur 1	Bur 2	Bur 3	Bur 4	Bur 5	Bur 6	Bur 7	
Straight	3.5		Pilot	VT1828 (2.3)	VT2838* (3.3)	—	—		Pilot	VT1525 (2.0)	VT1828 (2.3)	VT2838* (3.3)	—	—	—	—	
Straight	4.0		Pilot	VT1828 (2.3)	VT2838 (3.3)	VS3238* (3.5)	—		Pilot	VT1828 (2.3)	VT2838 (3.3)	VS3238* (3.5)	—	—	—	—	

*Denotes implant placement.

*Clinician experience and judgment should be used in conjunction with the Densifying Reference Guide recommendation and suggested use protocols.