




STOP Drilling Away Healthy Bone

Presenting Densah™ Bur Technology*: The Innovation That Makes Osseodensification Possible

Introducing Densah Bur Technology for implant osteotomy preparation from Versah™ LLC. Densah Burs have a non-excavating proprietary flute design that, when rotating at 800 – 1500 rpm in reverse, densifies bone. This technique, known as Osseodensification, autografts bone along the entire length of the osteotomy through a hydrodynamic process with the use of irrigation. When rotating clockwise, Densah Burs also precisely cut bone. The result is a consistently cylindrical and densified osteotomy. Consistent osteotomies and densification are important to implant primary stability and to early loading.

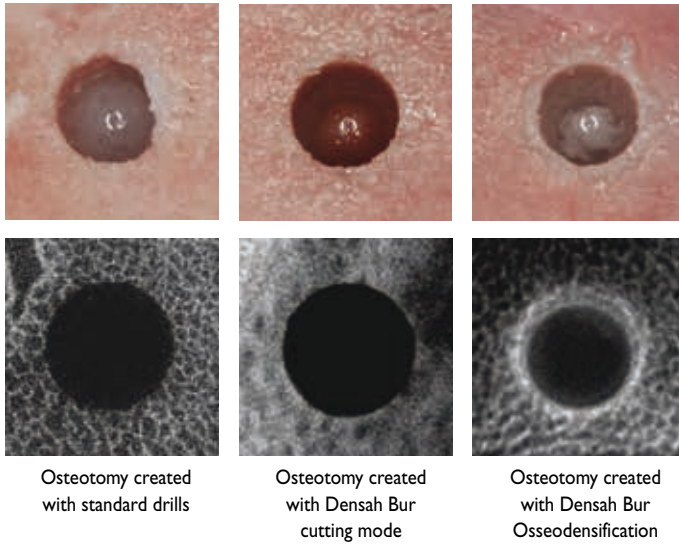

To order the newest innovation in implant dentistry, contact a Versah Customer Service Professional at 844-711-5585 or visit www.versah.com

